

Recoverable metal microdrive


Design name: metal_v7

Travel distance: 5.7 mm


Shell base: 3.1 x 5 mm (WxL)

More information is needed:


- contact me (Misi Voroslakos) directly at voroslakos@gmail.com

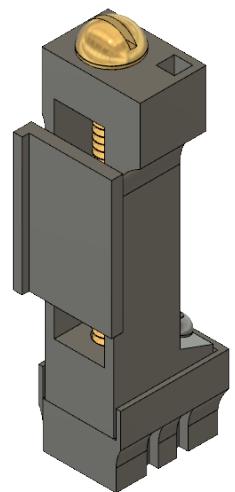


arm/shuttle

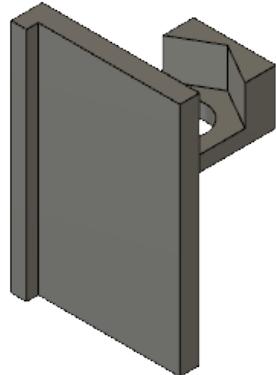

drive body

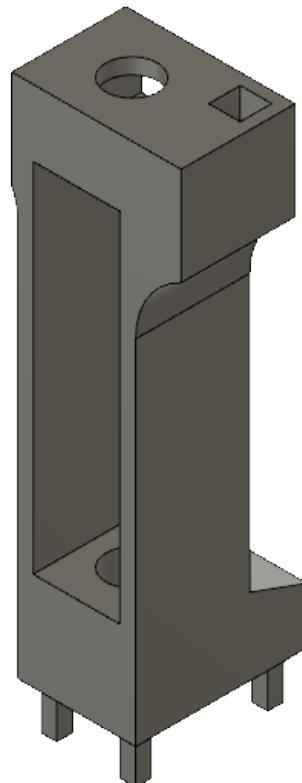

base

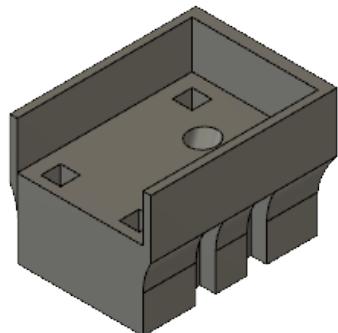

00-90 screw
1/2"



brass hex nut




000-120 screw
1/8"




male header pin

1. Tap base (000-120 tap).

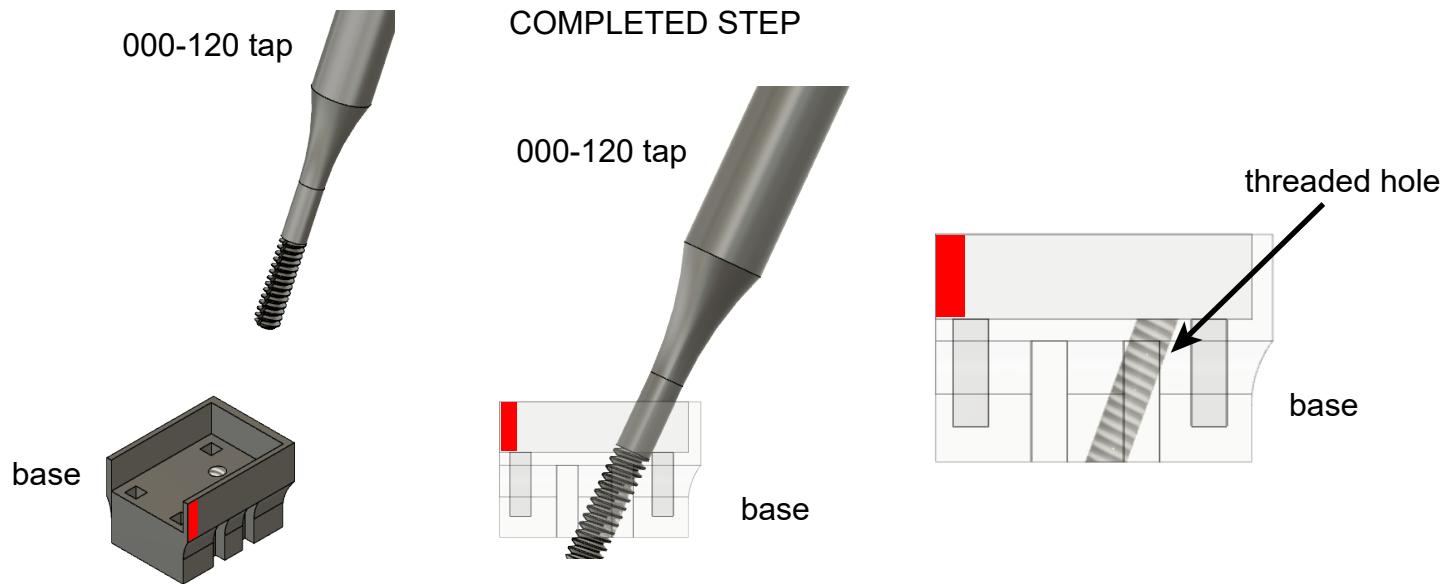


2. Insert and glue nut into arm.

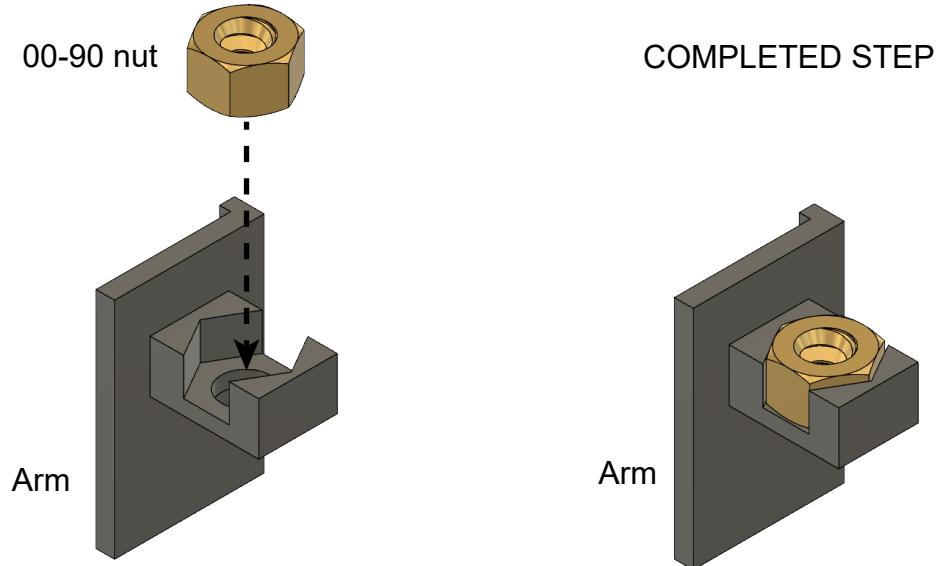


3. Solder 00-90 hex nut to 00-90 screw.


3a - Insert nut into drive.

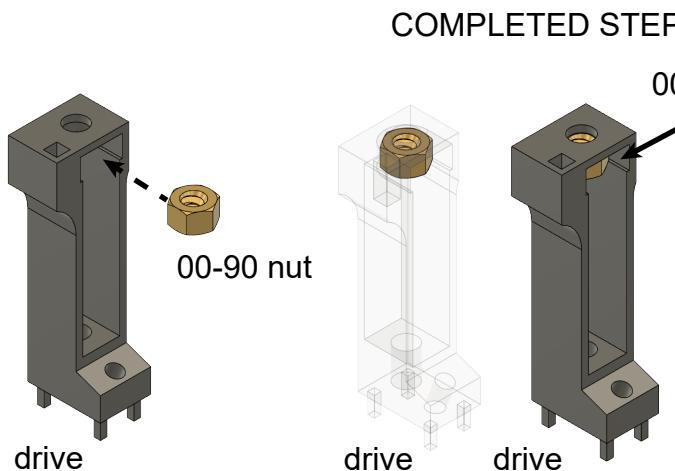


3b - Insert arm into drive.

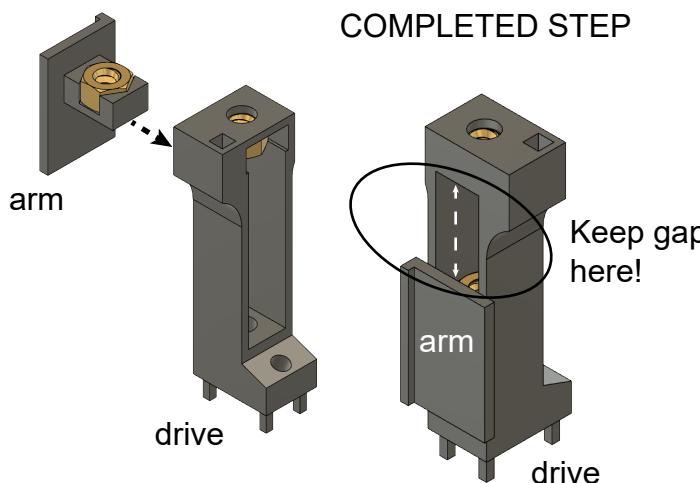


3c - Insert 00-90 screw into drive and tighten nut.

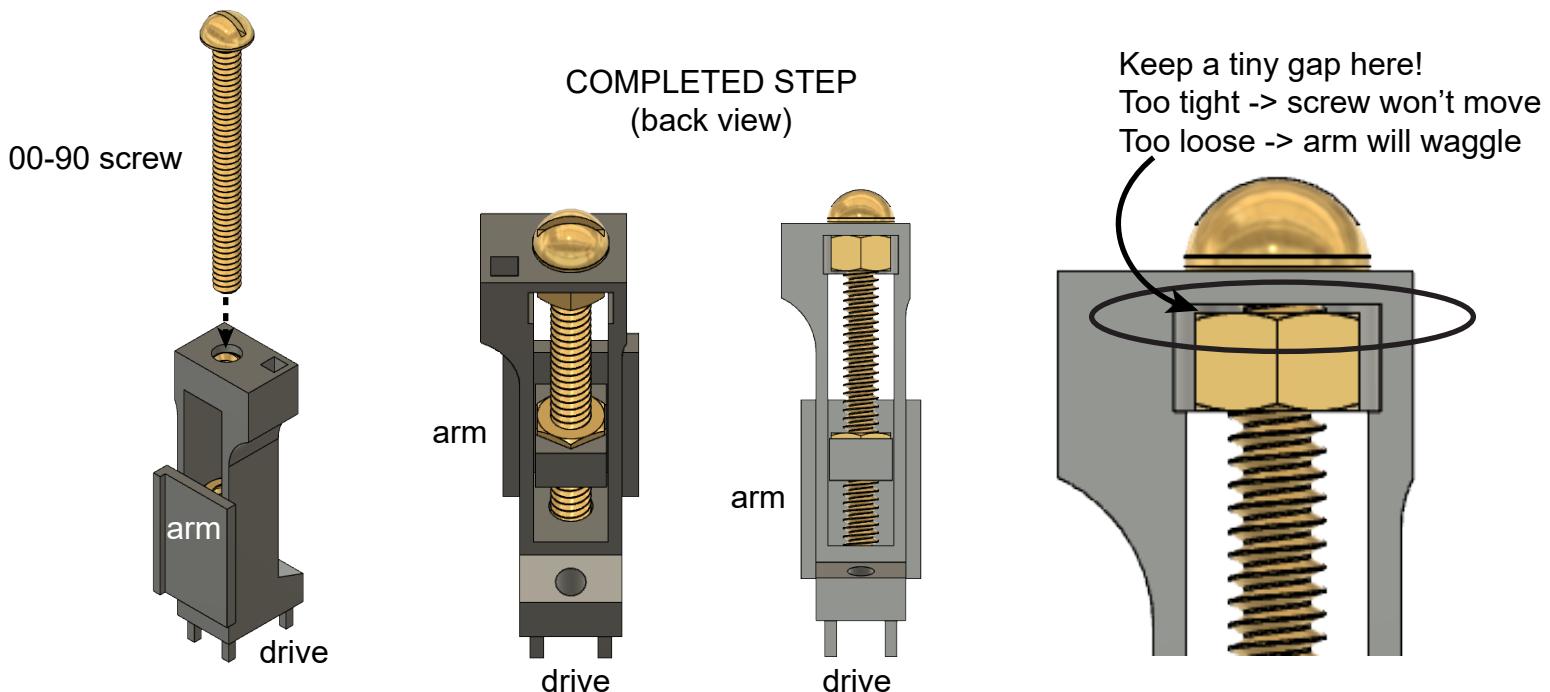


3d - Solder 00-90 nut to 00-90 screw.

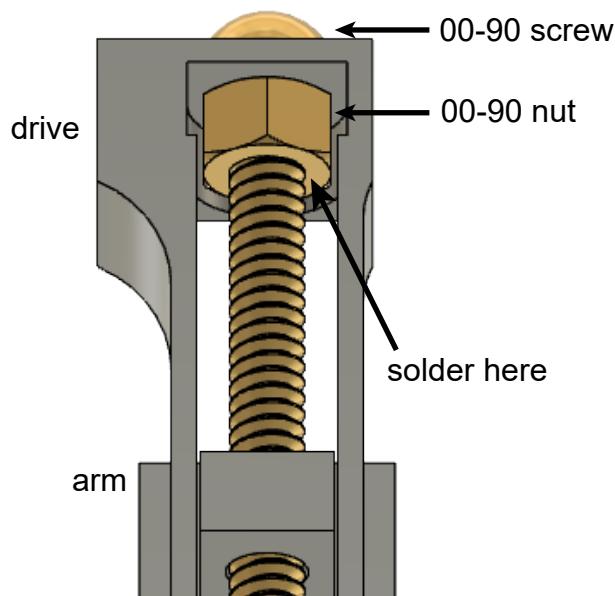


COMPLETED STEP

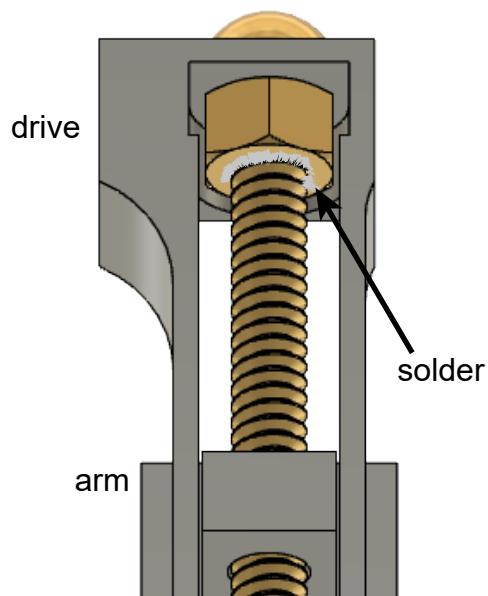

4. Attach base to drive.

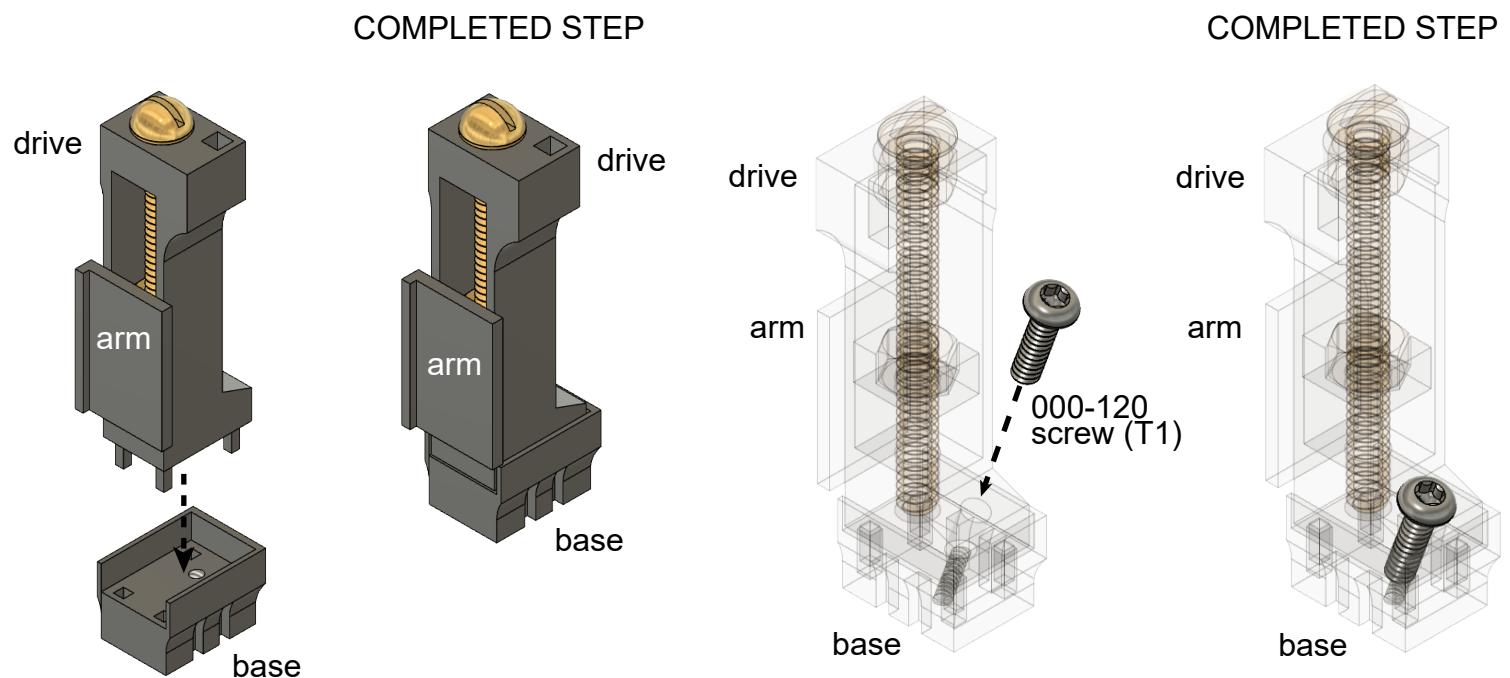

5. Insert male header pin into drive (if used during surgery).

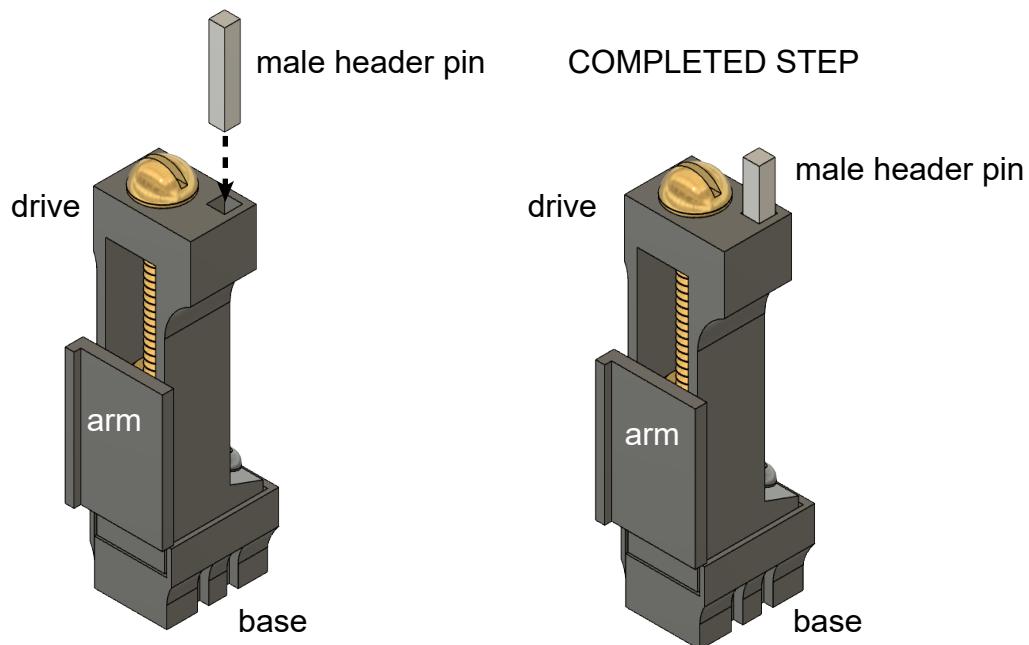


Implantation tool for recoverable metal microdrive


Design name: stereotax_attachment_metal_v7
drive_holder_metal_v7

More information is needed:

- contact me (Misi Voroslakos) directly at voroslakos@gmail.com

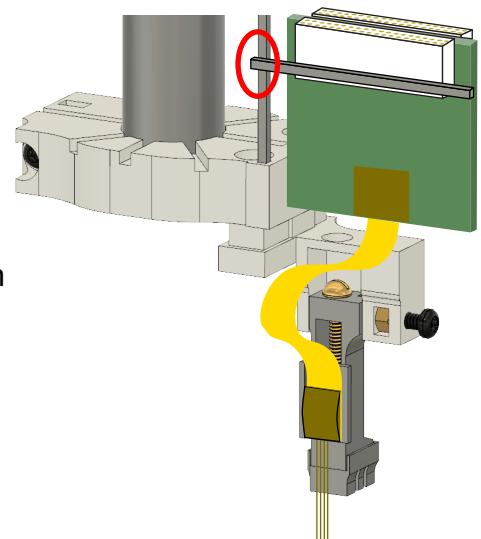

stereotax attachment

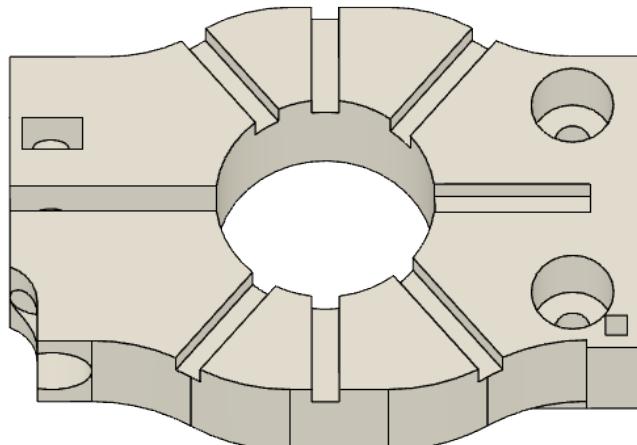


male header pin



00-90 nut
(4x)




drive holder

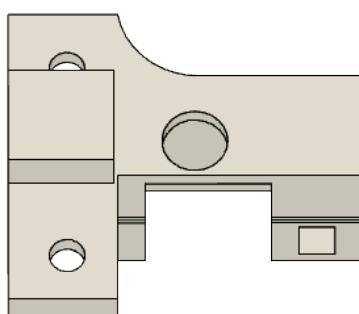


00-90 screw
1/4" (3x)
3/16" (1x)




1. Insert and glue 00-90 nut into stereotax attachment.

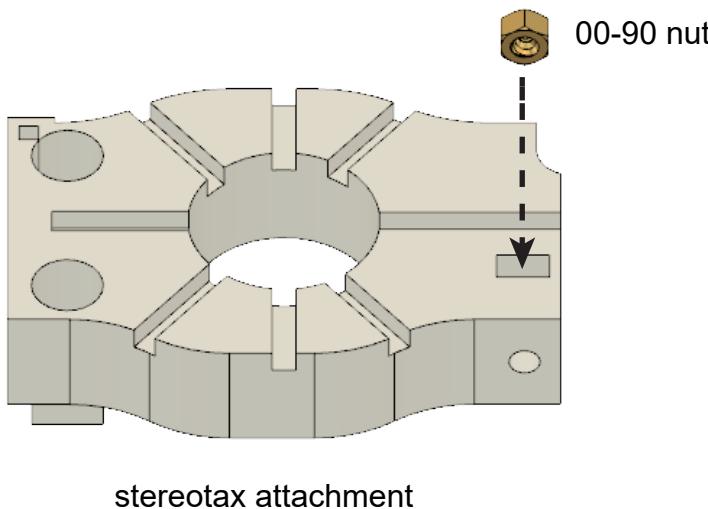


2. Insert and glue male header pin into stereotax attachment.

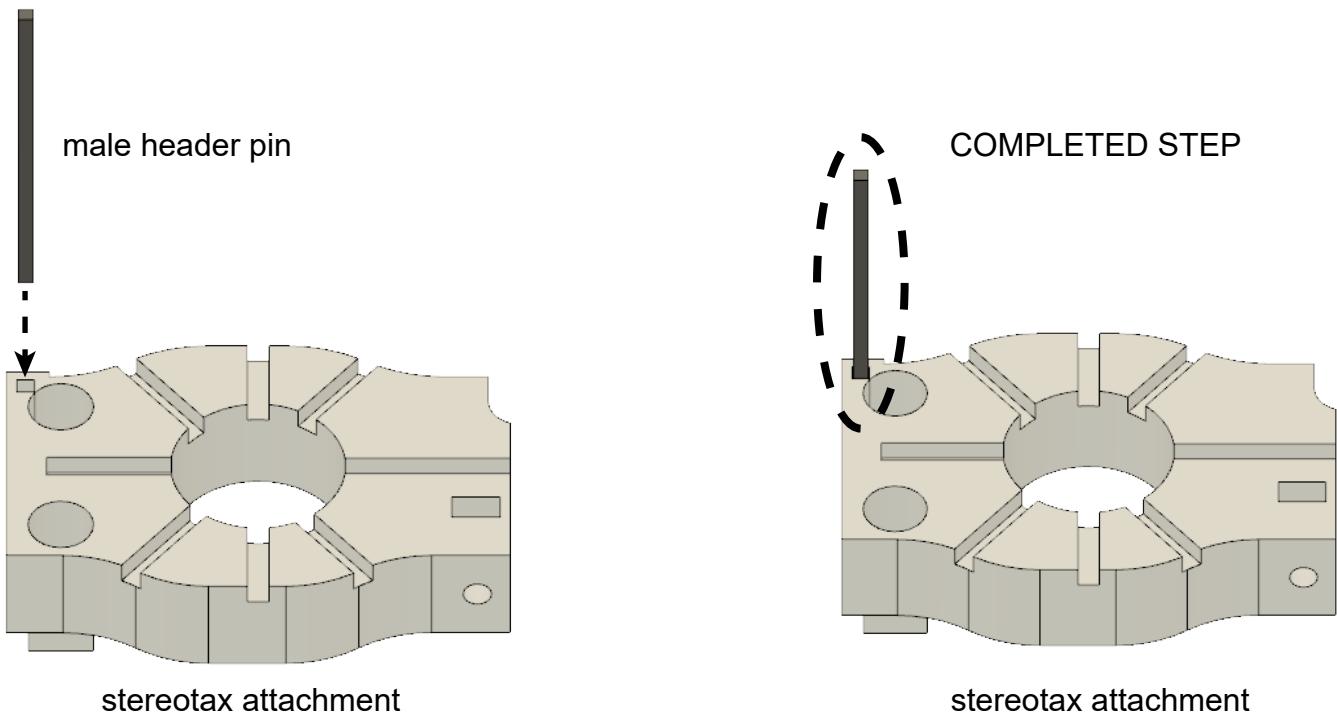


3. Insert 00-90 1/4" T2-screw.

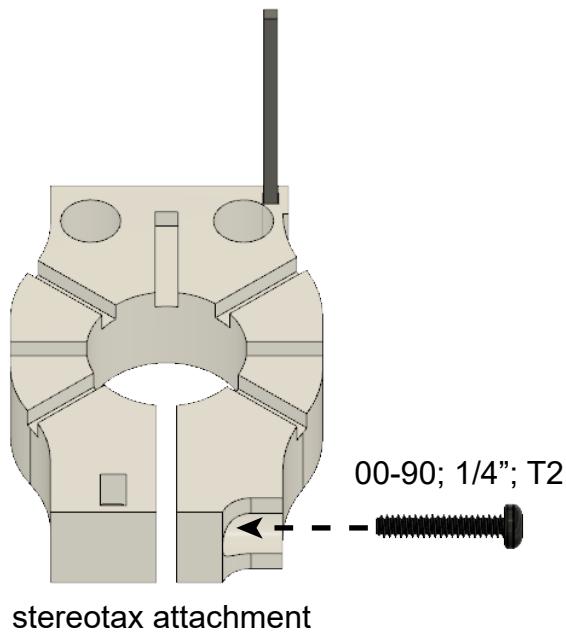


COMPLETED STEP

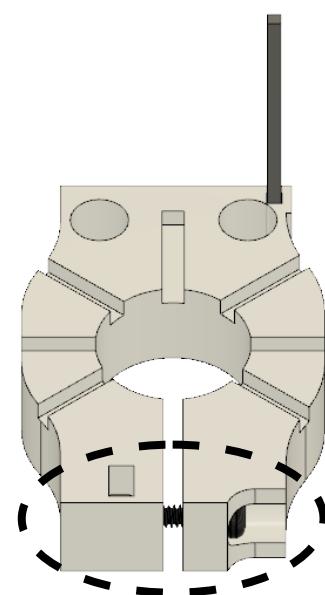

4. Insert and glue 00-90 nuts (3x) into drive holder.

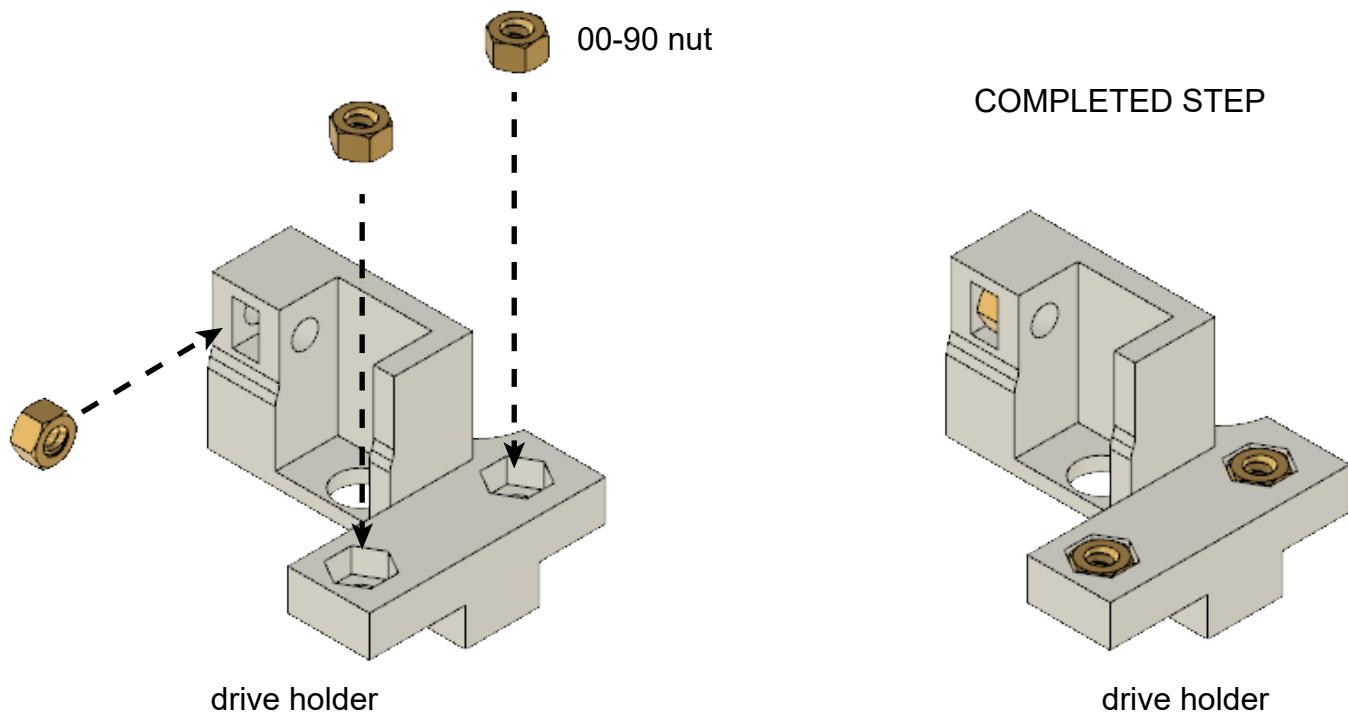

5. Insert 00-90, 1/4" T2-screws (2x).

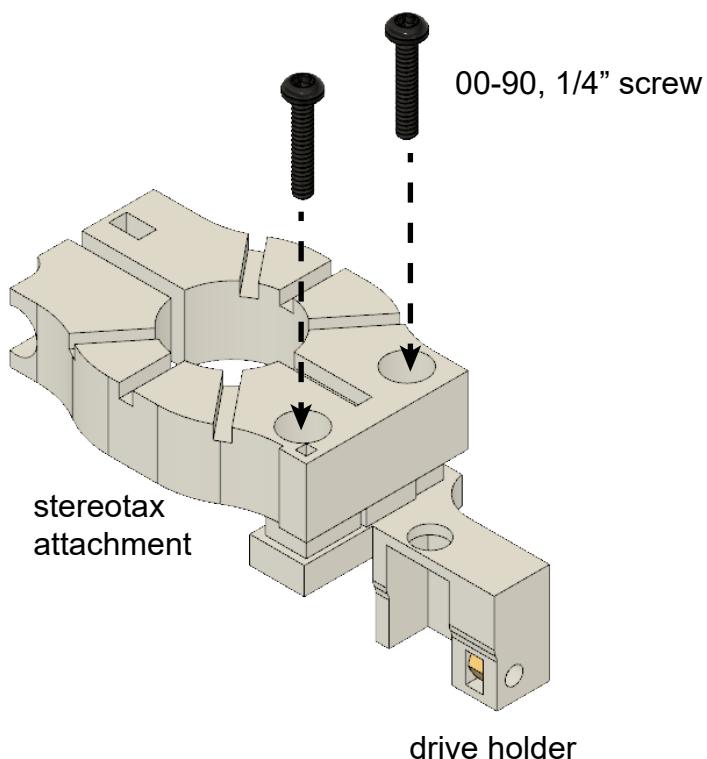

6. Insert 00-90, 1/4" T2-screw.

